

U.S. FDA

U.S. Food and Drug Administration (FDA)
approved medical device (510K 072534)

J.ELECTRIC SAFETY

Permitted by the Japanese Electric Safety Authority
(09-09087182)



Good Manufacturing Practice (GMP)
approved manufacturer (No. 05-794)

ISO 13485

Approved as an international
medical device manufacturer (08 0037 S.J)

ISO 9001

Approved as an international
electric equipment manufacturer (SQN1454/94740)

CE

International Electrical
Safety Approval (Eui-si-2009-0555)

KENCO

Korean Electrical Safety Corporation (KESCO)
Electric Appliances Safety Certificate (H407291-9001A)

KISTI

Electromagnetic Shielding
Safety Approval (LP2007-0258)

BBB

Credit Rating (RICHWAY2011)



Richway biomat 7000mx

Advanced Safety System

Natural Amethyst Infrared rays and Negative Ions



All people deserve to be cured of their wounds,
rejuvenated from the aging process, spared from illness
and saved from suffering. It is our Birthright!

Richway International Inc. The name you can trust!



Richway International Inc.
The Life Corporation
Richway International Japan, Co.
Fuji Bio Science Inc.
Richway & Life Corporation



A STATE OF THE ART

Amethyst BioMat Professional

It is designed to bring health and happiness to people.



RICHWAY



This booklet is the user guide for the BioMat. It describes in detail how to operate the BioMat, what to be careful with, and how to use the BioMat to get the best results.



THE
RICHWAY
GROUP

RICHWAY International Inc. 1314 South King street 520 Honolulu HI 96814 USA. Tel: 808-589-2800 www.richwayusa.com



R I C H W A Y


LIFE WITH A BIOMAT

History of Heat Therapy

The history of raising body temperature for therapeutic purposes can be traced back twenty four hundred years ago.


“Give me the power to produce fever, and I will cure all disease.”

Hippocrates (BC 460)

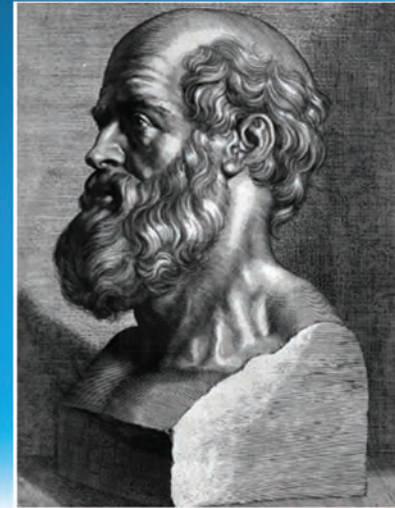



If I possessed the power to generate heat, I could cure every disease in the world. (Hippocrates, 450 BC)

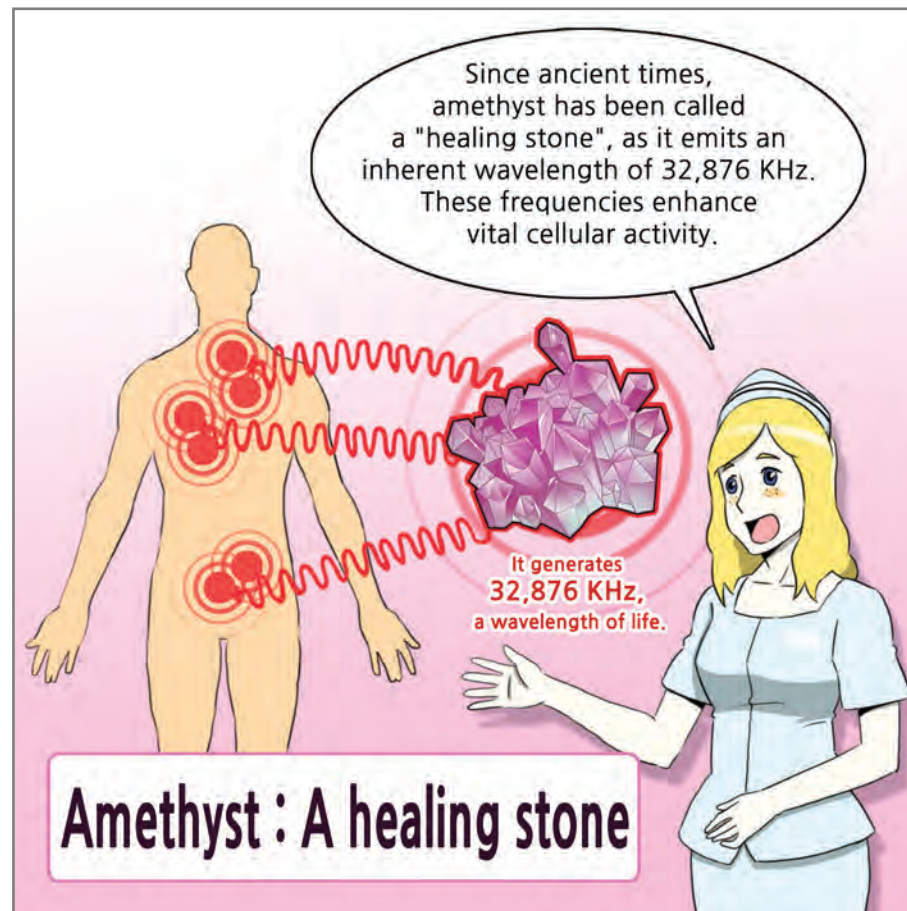


Far-infrared radiation is the safest energy in the world, as well as the light of life. (Hwang Sang, MD, 1995)

Type of light based on wavelength

All light has its own inherent wavelength. Gamma-rays being the shortest wavelength, with wavelengths increasing in sequence : X-rays, ultraviolet rays, visible rays (0.4~0.7nm) and then finally far infrared rays with deep impact and powerful thermotherapeutic energy.

"Amethyst emits more powerful far infrared rays with deep impact power and more powerful thermotherapeutic energy."

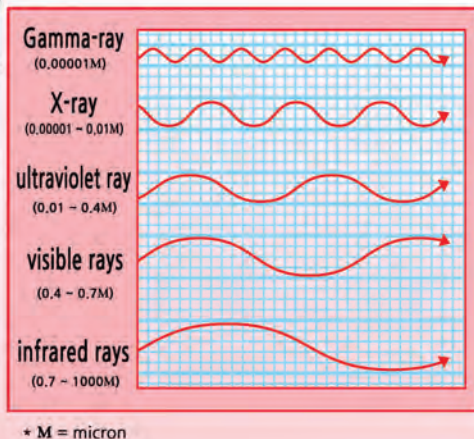

Heat

Far infrared rays from amethyst

What is the difference between heat and far infrared rays?
It is difficult for external heat to penetrate into our body, but far infrared rays can penetrate deeper into our bodies and generates intense frictional heat which stimulates our cells to move.

Photos showing the difference between an ordinary mat and the BioMat with far infrared rays.

When lying on an ordinary mat

When lying on the BioMat

Easy, isn't it?

A good surgeon has an eagle's eye, a lion's heart, and a lady's hand. (British Saying)

Differences between ordinary far infrared rays and far infrared rays from amethyst


far infrared rays with long wavelengths (which are disruptive and cause emotional disorder)

far infrared rays with shorter wavelengths (6.5 - 12 micron) (which have deep impact power and are quite comfortable)

"6.5 to 9.32 micron far-infrared radiation with the human resonance wave"

Active movement of cells generates frictional heat and raises your body temperature. Especially the amethyst far infrared rays of the BioMat which have a shorter wavelength (6.5 - 12 nm). They have deeper impact power and generate heat more powerfully and comfortably than that of ordinary far infrared rays.

The BioMat has a Hideo figure design.

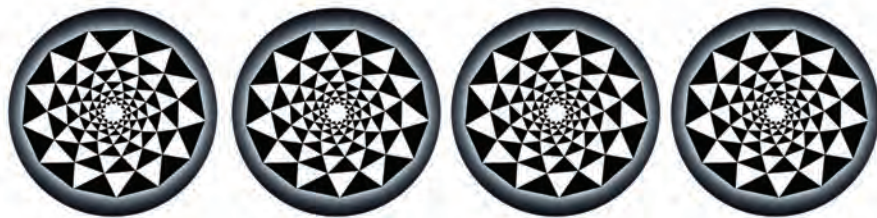

The Dewa Hideo figure is famous for its bio-resonance energy emission.

"6.5 to 9.32 micron far-infrared radiation with the human resonance wave, which is very effective in treatment"



The BioMat is designed to maximize bio-resonance energy as a result of its Hideo figure.

Exercise is enjoyment. There is no happiness without working out. (ArnoldSchwarzenegger)

Excellent heat insulation effects of the BioMat

The BioMat warms up your body and keeps it warm as the far infrared rays change your physical condition generating heat and stimulating the mitochondria in your cells.

"6.5 to 9.32 micron far-infrared radiation with the human resonance wave, which has excellent heat retention properties"

The BioMat test data of Telesum in the US

Set low temperature in back of the head (113°F/45°C)

Set low temperature in middle part of the waist (113°F/45°C)

time	temperature	time	temperature
0 second	24.4°C	0 second	24.4°C
30 second	26.6°C	30 second	26.6°C
60 second	24.5°C	60 second	24.5°C
90 second	26.6°C	90 second	26.6°C
1 hour	27.6°C	1 hour	27.6°C
Check change of body temperature after turning off the BioMat.			
2 hour	26.1°C	2 hour	33.9°C
4 hour	26.7°C	4 hour	33.8°C
6 hour	27.9°C	6 hour	34.5°C
8 hour	27.3°C	8 hour	34.7°C

*Body temperature has maintained the same level 8 hours after turning off the BioMat.

Isn't it amazing?

Before using the BioMat.

Two weeks after using the BioMat.



Above is a photo of a patient with a blood circulation problem. After using the BioMat, you can see that his blood circulation has been improved considerably. It is because the patient has changed his physical condition to maintain high body temperature by himself.

"6.5 to 9.32 micron far-infrared radiation with the human resonance wave, which promotes the circulation of the blood."

Dr. George Grant

There are about 1000 times more toxic heavy metals used than 50 years ago in the system of modern men. (The Journal of the American Medical Association, 1974)

I really hate far infrared rays ~

I want to live in a world without the BioMat ~

Me too, me too!

Harmful substances

Amethyst far infrared rays will change your physical condition to raise your body temperature.

One of the prominent phenomena of hyperthermia treatment with the BioMat is the change of the physical condition to keep high body temperature by itself. Such change of the physical condition facilitates the blood circulation and activates more than 3000 essential enzymes, so harmful materials in the body will be discharged naturally.

Blood purification from amethyst far infrared rays

1

Before using the BioMat

2


30 minutes after using the BioMat

This photo of blood was taken after using the BioMat for 30 minutes. You can see that turbid blood with coagulated red blood cells was changed to healthy and ideal blood suitable to carry oxygen thanks to the far infrared rays of the BioMat. Isn't it amazing?


The average adult battles cancer furiously about six times in his or her lifetime.

The BioMat hyperthermia function will enhance cardiovascular activities

Before using the BioMat(2010.5.29)

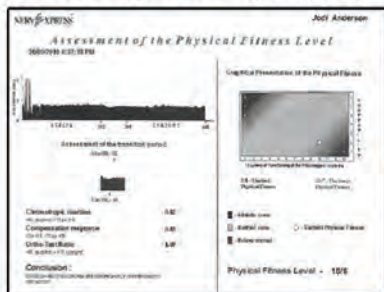


30 minutes after using the BioMat(2010.5.30)

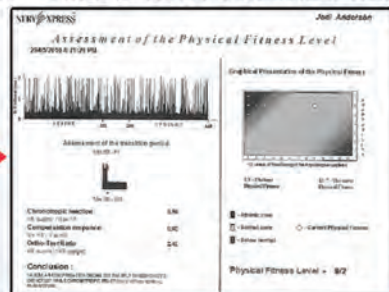

This photo shows the heart of a 57 year-old woman. Her health condition was very poor with reduced cardiovascular functionalities and arrhythmia. However, only three minutes after using the BioMat, her cardiovascular system was greatly enhanced. Likewise, we can confirm that the hyperthermia function of far infrared rays will affect our circulatory organs considerably.

Your immune system will deteriorate if your body temperature drops.


All diseases come from reduced body temperature. Such lowering of body temperature is called "Cold hypersensitivity". A symptom of cold hypersensitivity is getting cold hands, feet, & other areas that are far from the heart, and then an overall drop in body temperature.

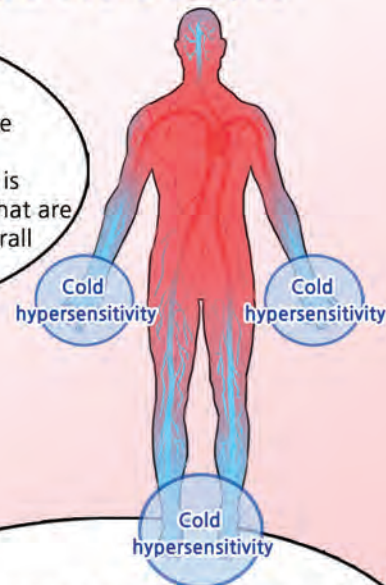

Especially, when the temperature of hands and feet drops due to the lowering of outside temperature, terminal capillary vessels are contracted, so overall blood circulation will be hindered.

"Kill off the bad bacteria and clean up the waste"; that's what your immune system works to do all day.


If your body temperature drops by 1°C, your vital metabolic activities will be reduced.

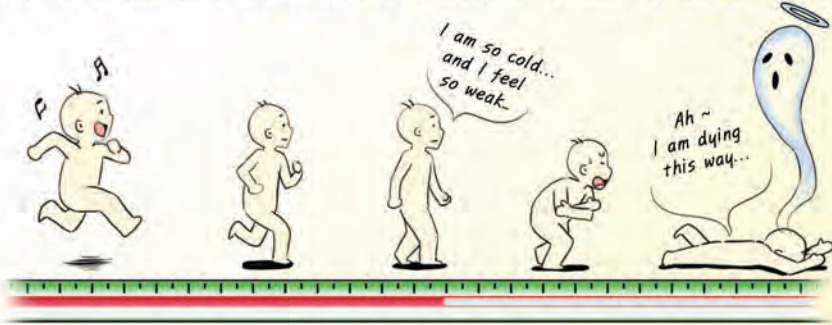

Even if your body temperature drops by only 1°C, more than 3000 essential enzymes cannot perform their function for your body.




If your body temperature rises by 1°C, then enzyme activities will be increased about 40%

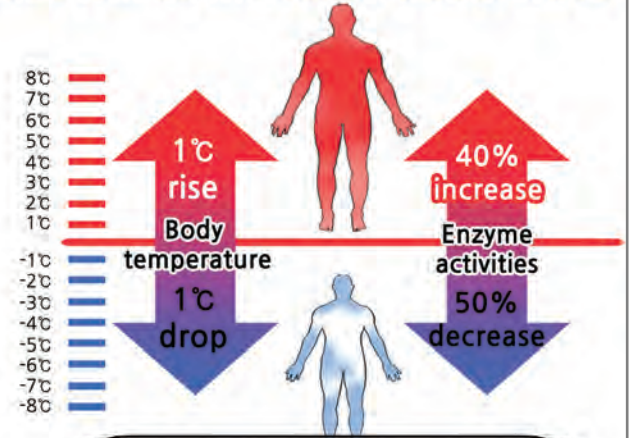

If your body temperature rises by 1°C, then enzyme activities will be increased by more than 40%. On the other hand, if your body temperature drops by 1°C, then enzyme activities will be decreased by more than 50%. In other words, people with cold hypersensitivity tend to get sick easily due to weakened immune systems.

The immune system seeks out and destroys all the cells that don't assist in performing normal physiological activities.


Your body temperature shall be over 36°C on average.

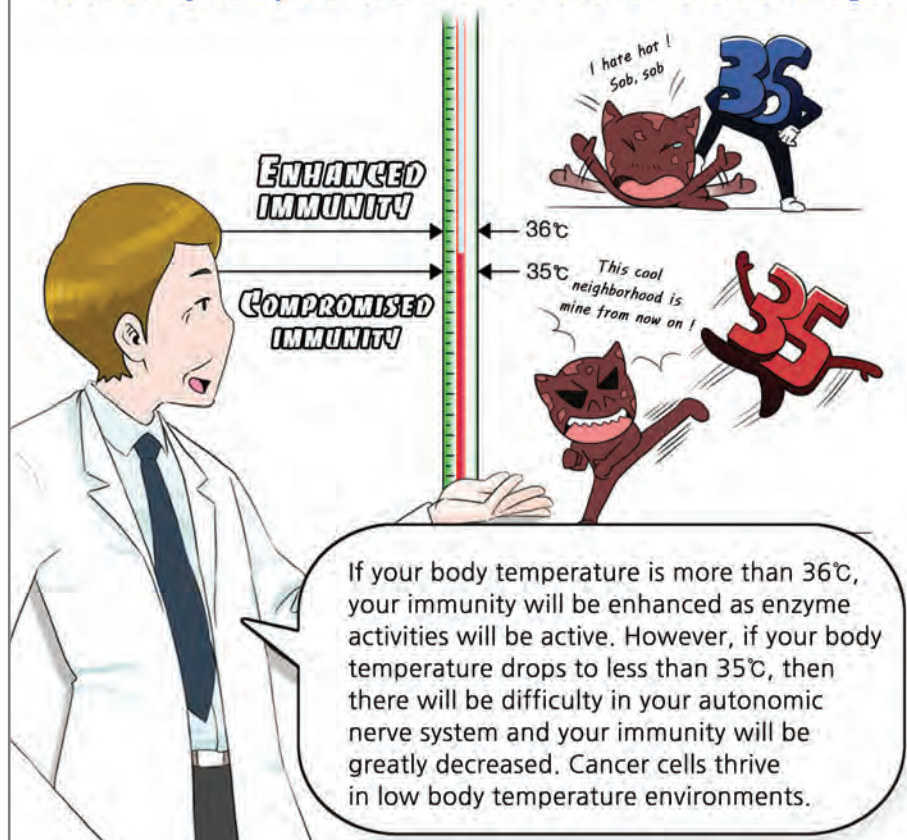


Report of Harvard Medical School



Enhanced enzyme activity due to an increase in body temperature


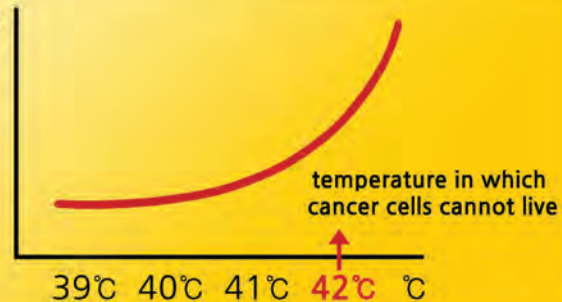


Harvard Medical School found out that cancer cells are vulnerable to heat, and will be destroyed at temperatures over 42°C. The rise in body temperature will weaken activities of cancer cells.

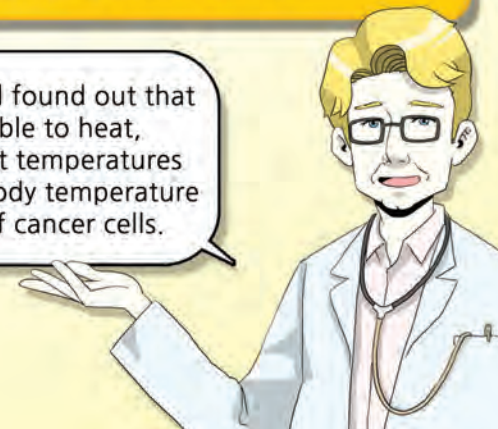


Every day your cells are divided about two billion times.

Life with a Biomat.

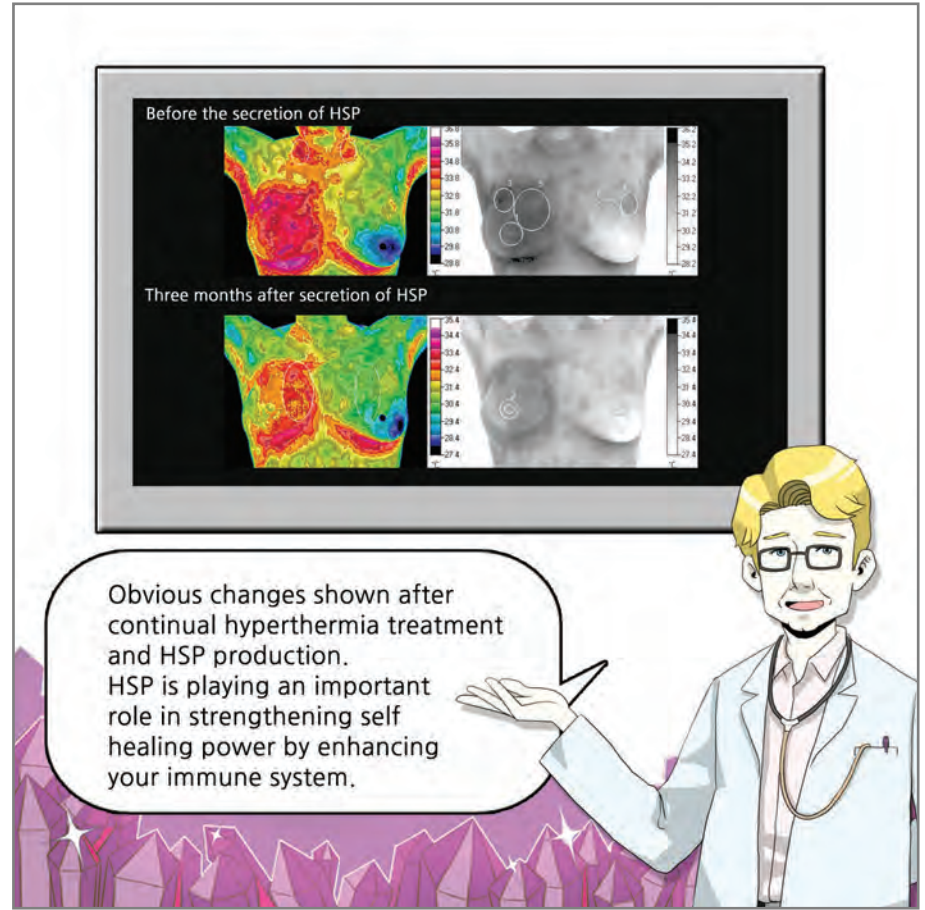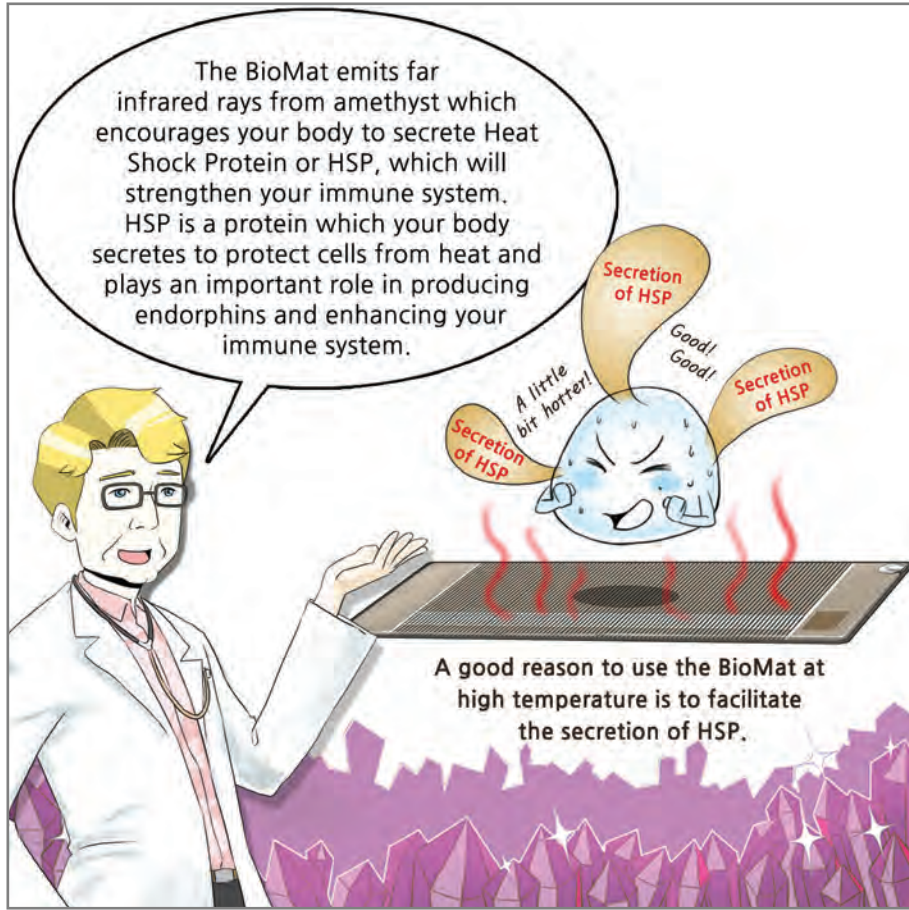

Eating a diet that constantly raises blood glucose levels is like trying to put out a forest fire while someone nearby is throwing gasoline on the trees. (Lorraine Day, MD, 1984)

HSP increases the production of endorphins, NK-cells, T-cells, and lymphocytes.

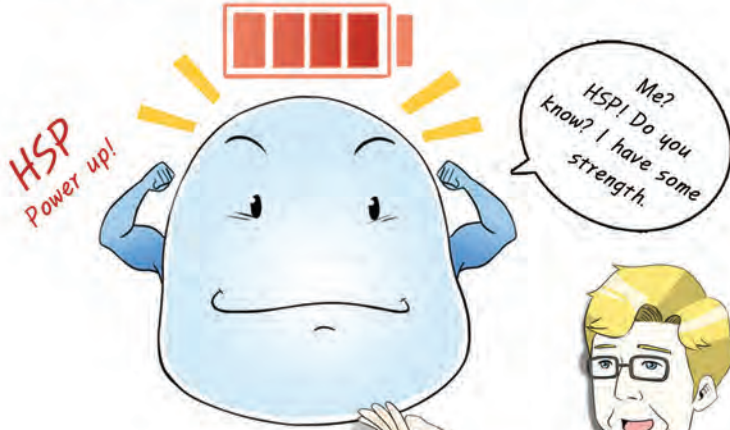

HSP
Power up!



Me?
HSP! Do you know?
I have some strength.

1. Enhances Immune System.
2. Reduces Fatigue.
3. Relieves Pain.

The BioMat reduces stress hormones by 78%.

Magnetic resonance imaging test, Electroencephalogram test, Electrocardiogram test, Blood cortisol test


Level [mcg/dcl]	Pre Test Results Before use the Biomat	Post Test Results After use the Biomat
Magnetic resonance imaging test [0-100]	68	34
Electroencephalogram test [100-1000]	575	285
Electrocardiogram test [1-10]	2	6
Blood Cortisol test [5-25]	16	7

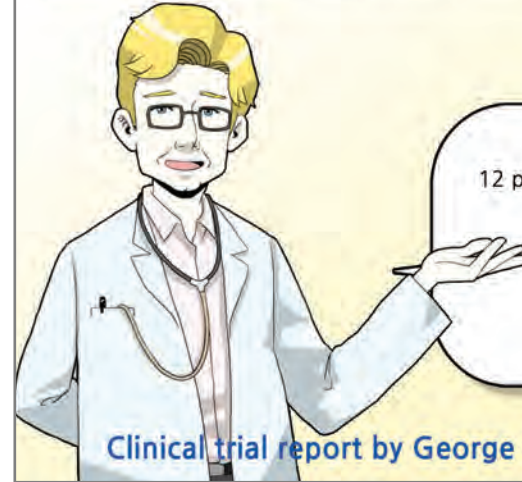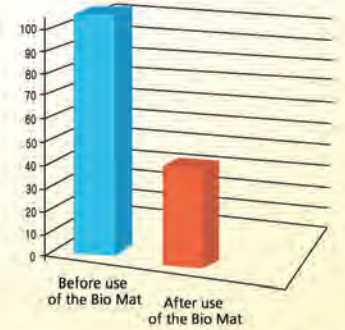


After an experiment in which 12 patients used the BioMat for one hour a day to find out if the BioMat can reduce stress levels efficiently, the results showed that stress hormones, including cortisol, are reduced by 78%.

Clinical trial report by George Grant, a Canadian M.D.



If you can't be free from stress, it is like fighting a disease with your hands tied.

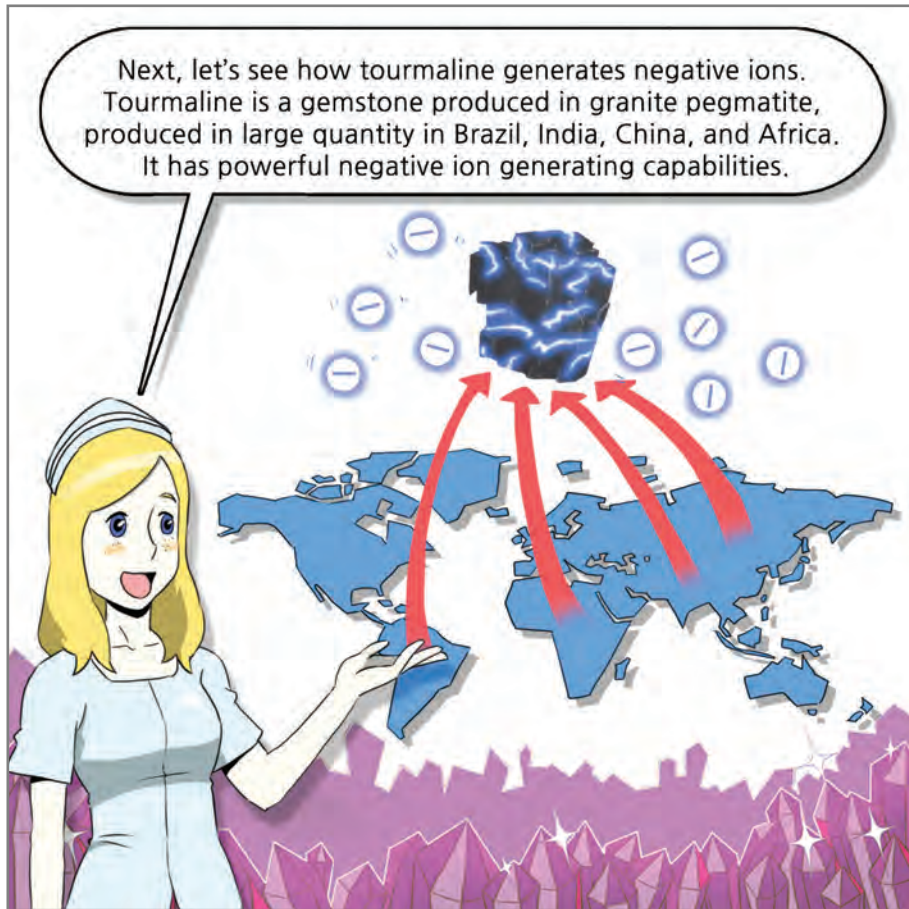

People put stress upon themselves. One may as well altogether dispense with feelings including hatred, abomination, obsession, and selfishness.

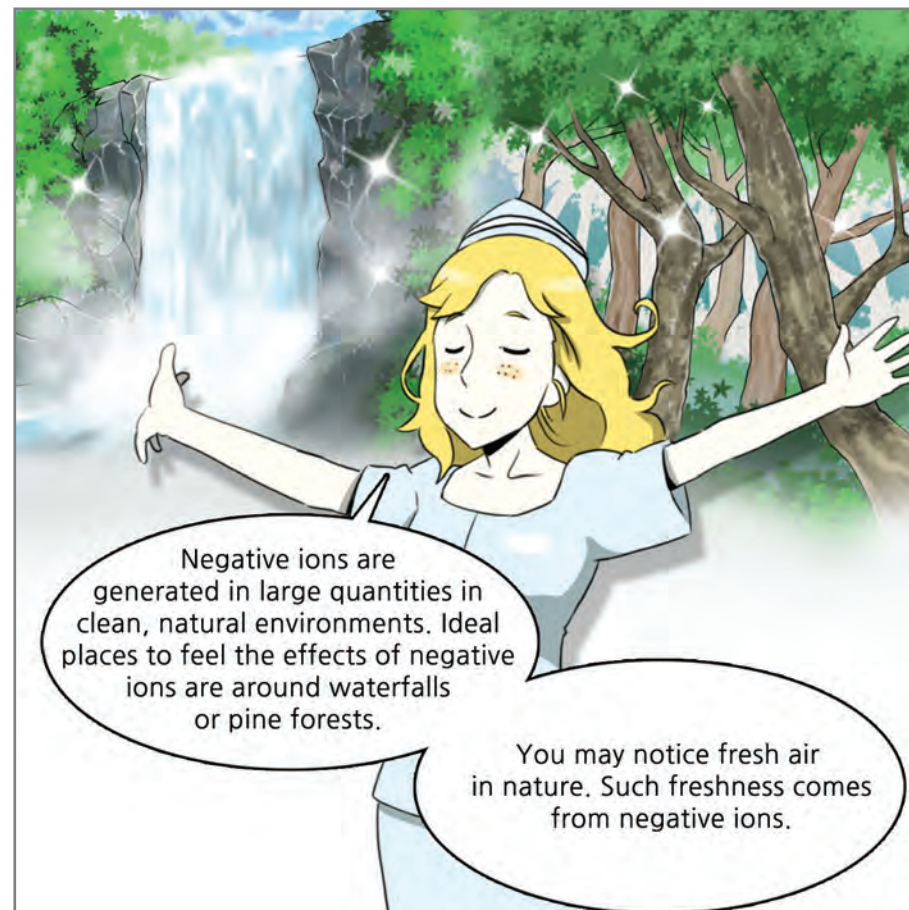

Only after you have embraced the feelings of delight, passion, gratitude, and consideration, will you receive the life from God.

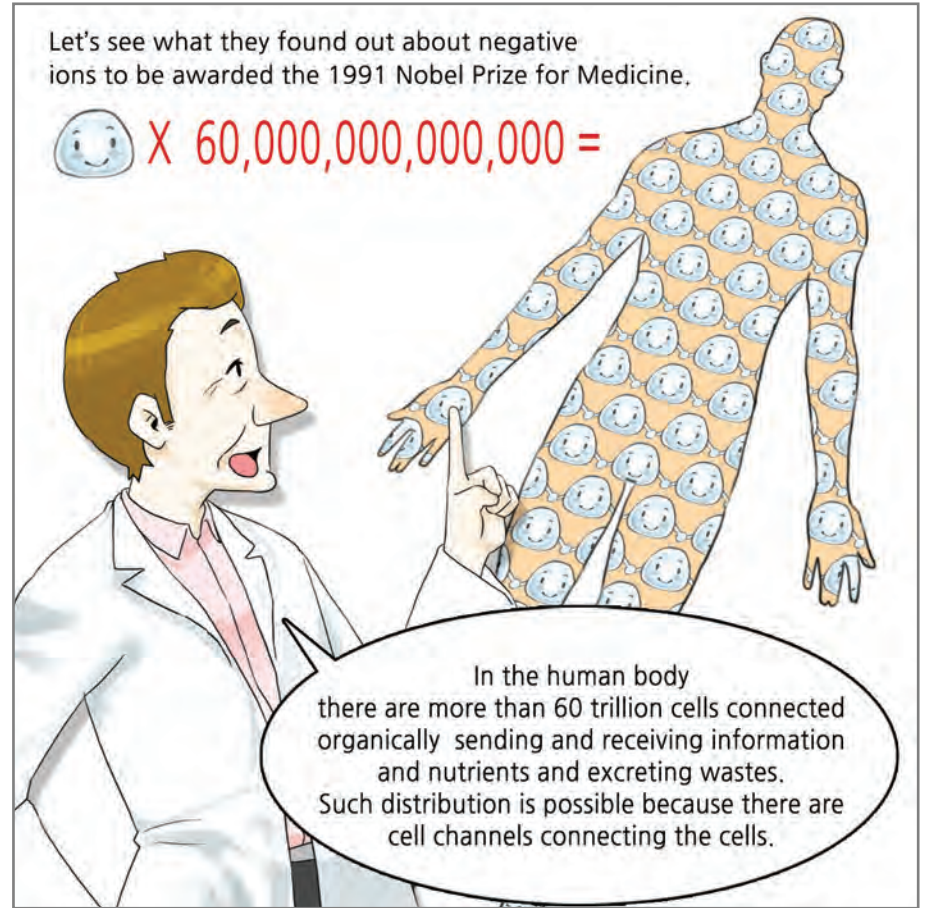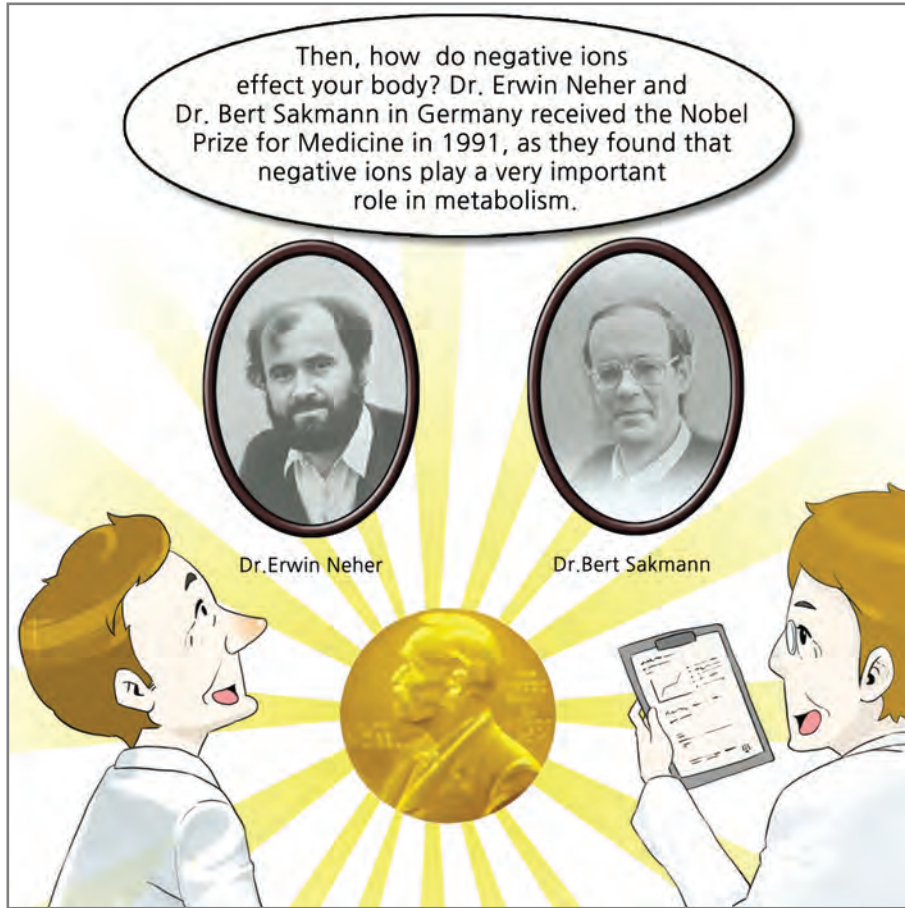

Enjoy beauty, music, laughter, and enjoyment around you.

Without negative ions, your cell channels do not open.




However, if there are no negative ions in cells, then cell channels will be closed. Thus, nutrients cannot be passed to each and every part of your body, so activities of cells will be reduced creating pain.

Wastes that aren't excreted accumulate in the body causing adverse health.



You feel pain when your brain detects electrical signals of positive ions. Negative ions generated by the BioMat will be absorbed into your body through your skin and your breath. They will facilitate cellular metabolism, enhance your vitality, reduce pain and distribute various nutrients and energy throughout your body. They will also cleanse blood and have a tranquilizing effect. Negative ions are also good for recovering from fatigue.

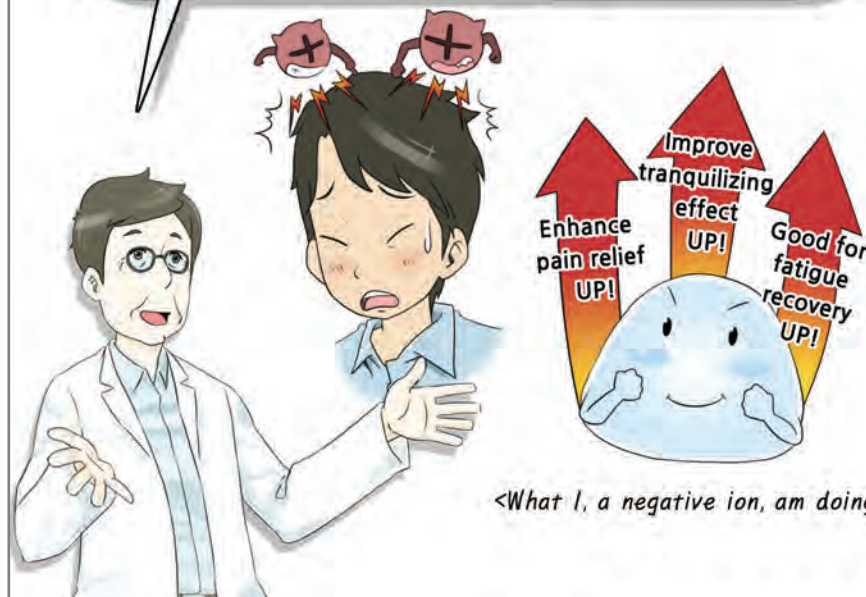

<What I, a negative ion, am doing!>

In the end, the love you take is equal to the love you make. (Abbey Road, 1969)

Components of the BioMat

< Mat >



< Controller >



< Quilted Pad >



< Pillow >



< Carrying case >



* Pillows sold separately.

Now,
let's see what's in the box.
A mat, a controller, a pad, and
a carrying case.
You can also purchase a pillow
separately.

Types of BioMats

King and queen



Single



Professional (BMP)



Mini



The BioMat is available in
king, queen, single,
professional (BMP), and mini.
King, queen and single are
designed for standard beds,


while professional (BMP) size is made
for intensive hyperthermia treatment and
is ideal for hospital or massage beds.
The Mini-sized mat can easily be carried
during travel or used to cover a chair.
It is also used to cover a patient during
hyperthermia treatment.

How to use the BioMat controller

Now~
let's see how to
use the controller.



A temperature indicator is shown on the display of the controller. Left side of the display shows a Celsius scale and right side a Fahrenheit scale. In addition, there's also a far infrared ray indicator, an electromagnetic interception indicator, and a built-in alarm system on the controller. On the right side, there's a power switch, temperature control buttons and time control buttons.

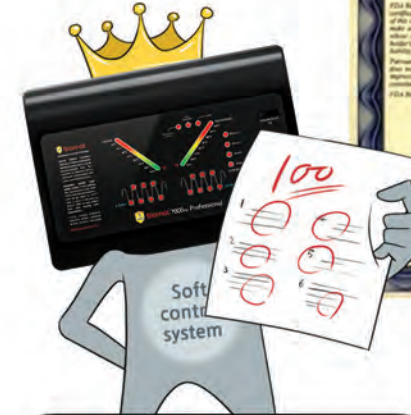

The BioMat controller is designed to safely control the generation of far infrared rays and negative ions with simple switch operation.



Diseases which cannot be treated through diet are incurable with any medicine. (Hippocrates, 450 BC)

First of all, connect the mat to the controller.




Connect the mini-mat to the controller.



Frustration and anger result in diseases. (French Saying)

With a king-size or queen-size mat,
you can control the temperature of the left
and right side independently.

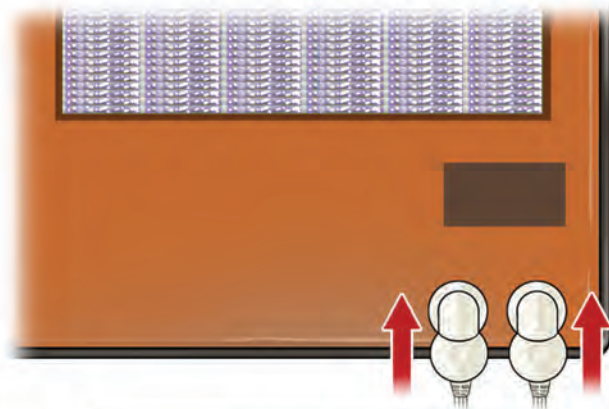


As there are two connectors on a king-size or a queen-size mat, you can control temperature of left and right side independently.

* You can control the temperature of the right and left side independently only when connecting two controllers, or with one double-sided king or queen controller.




Connection of the BioMat and the controller

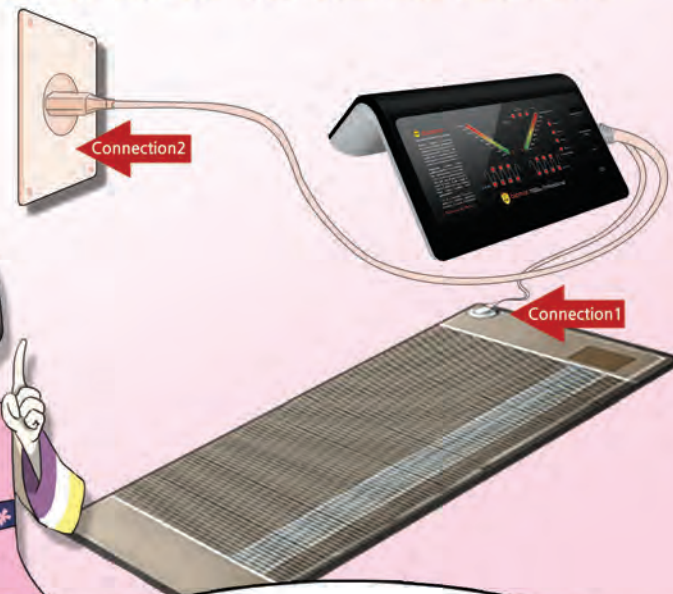

After connecting the BioMat and the controller, insert the electrical plug of the controller into an electrical outlet on the wall. Now you are ready to use the BioMat.


Give thanks in all circumstances. Gratitude is a medicine to treat your body and heart.

Always use the quilted cotton pad.

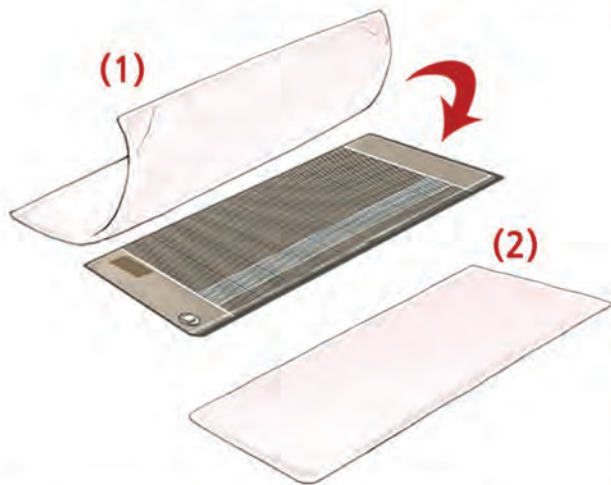


Whenever using the BioMat, always cover it with the supplied cotton pad. The pad has antibacterial properties and maximizes thermotherapeutic effects of far infrared rays.

**We recommend you cover your mat with the pad to maximize the hyperthermia effect all over your body.*

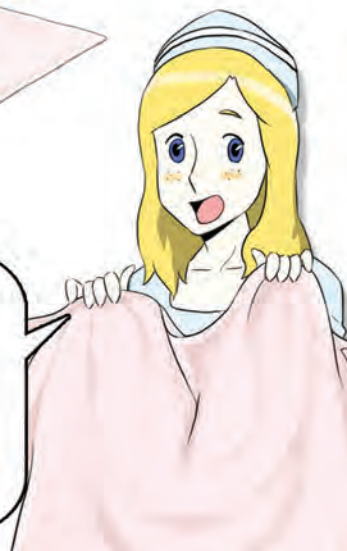

BioMat Power On 1

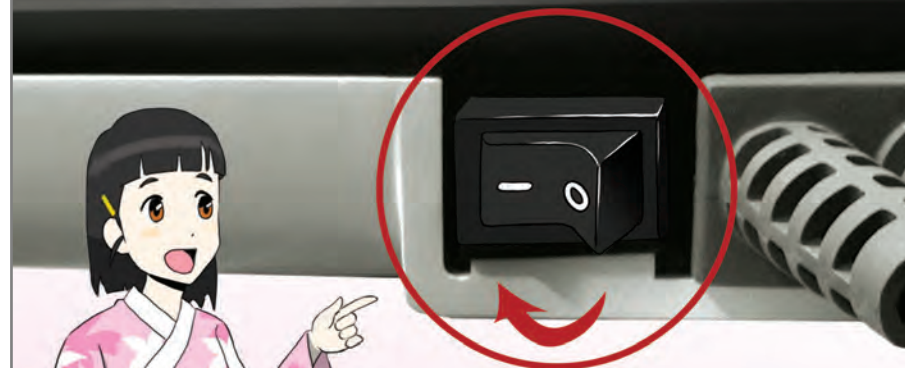

After connecting the BioMat and the controller, turn the power switch on the side of the controller to ON.

Forgive a person a day, and you can remove a disease a day. (British Saying)

How to control the generation of negative ions and far infrared rays

Temperature Indicator

Negative Ion Indicator

biomat 7000[™] Professional

If you push the power button on the right center console, only negative ions will be discharged without any other thermotherapeutic effects, and the negative ion indicator will be activated. If you push the temperature up button, the negative ion generation will be stopped and the far infrared ray thermotherapeutic function will be activated. When the mat reaches the selected temperature, it will discharge negative ions once again automatically.

How to set the BioMat timer


12 hours
8 hours
4 hours
2 hours
time setting

biomat 7000[™] Professional

Now it is time to set the timer. Default duration is 8 hours, however you can set 2 hours, 4 hours, 8 hours, or 12 hours.

Darkness is conquered not by fighting but by lighting.

How to set the BioMat temperature

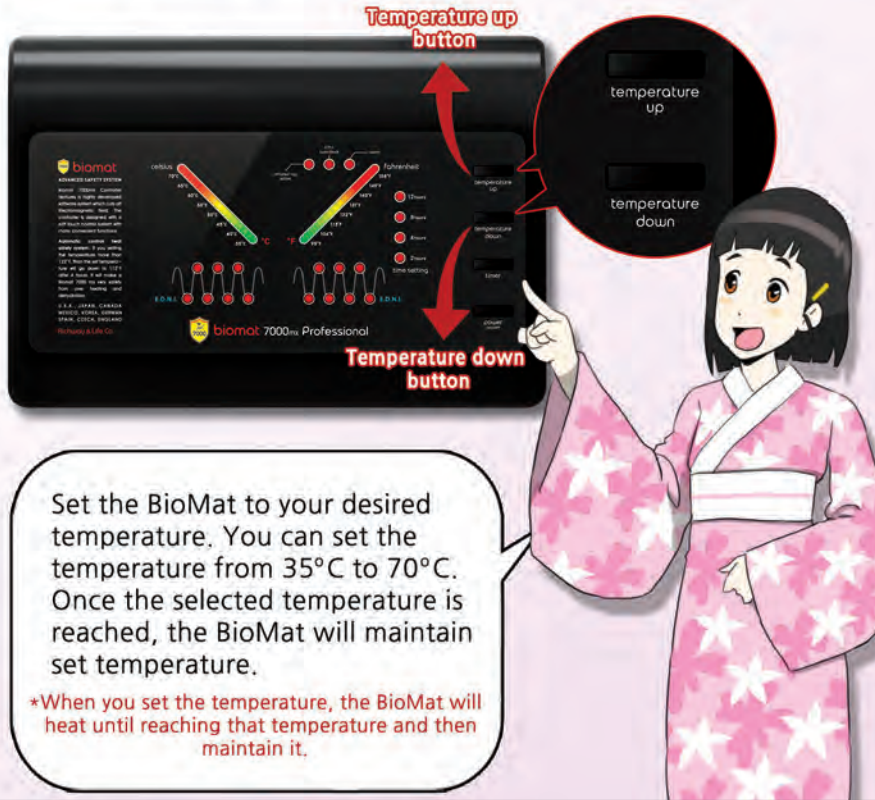

Temperature up button

Temperature down button

temperature up


temperature down

biomat 7000 Professional

Set the BioMat to your desired temperature. You can set the temperature from 35°C to 70°C. Once the selected temperature is reached, the BioMat will maintain set temperature.

*When you set the temperature, the BioMat will heat until reaching that temperature and then maintain it.

The BioMat's function of blocking electromagnetic waves

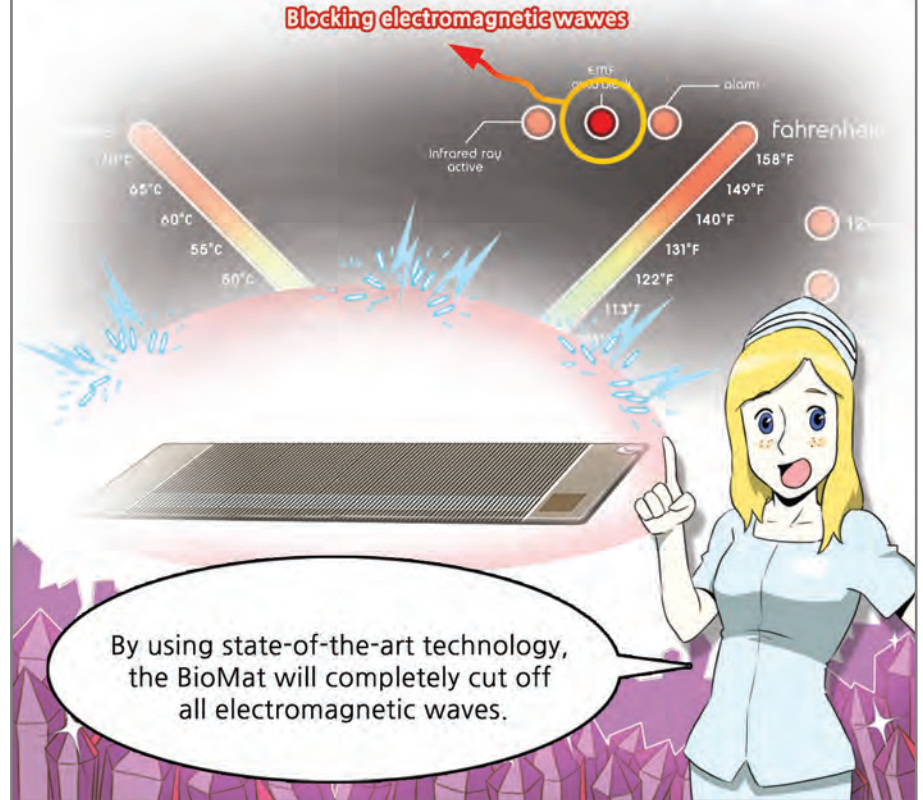

Blocking electromagnetic waves

Infrared ray active

alarm

Fahrenheit

158°F

149°F

140°F

131°F

122°F

113°F

111°C

85°C

80°C

55°C

80°C


By using state-of-the-art technology, the BioMat will completely cut off all electromagnetic waves.

Hyperthermia has been used since several thousand years ago, and doctors call it diaphoretic detoxification therapy.

Certificate of blocking electromagnetic waves

The BioMat, a medical thermotherapeutic device, has been certified by the KETI (an international organization for electromagnetic wave inspection) to have completely cut off electromagnetic waves.



K E T I Korea Electric Testing Institute	
EMC TEST REPORT	
Test Report No. : JP2007-0258	
Applicant Information	
Client	Rich Way & N.F
Client Address	ROOM6F, 66-4, BUN-DONG KANGBUK_GU, SEOUL, KOREA
Manufacturer	SU HYUN KIM
Manufacturer Address	ROOM6F, 66-4, BUN-DONG KANGBUK_GU, SEOUL, KOREA
EUT Information	
Equipment	Bio Mat professional Amethyst
Model/Type	BMP
Serial No.	N/A
Trade Mark(if any)	N/A
Rating	AC 120 V, 60 Hz
Receipt of the EUT	AUG. 28, 2007
Tested in period	AUG. 20, 2007
Date of issued	SEP. 4, 2007
Test Specifications	
IEC62233-2005 : Measurement methods for electromagnetic fields of household appliances and similar apparatus with regard to human exposure.	
Test Result	PASS
In the configuration tested the EUT complied with the standards specified above.	
Tested by (KETI in Korea)	J. E. Na / Tester
Approved by (KETI in Korea)	C. W. Park / Team Leader
<small>Remarks : This report details the results of the testing carried out on one sample, the results contained in this test report do not relate to other samples of the same product. The manufacturer should ensure that all products in series production are in conformity with the product sample detailed in this report.</small>	

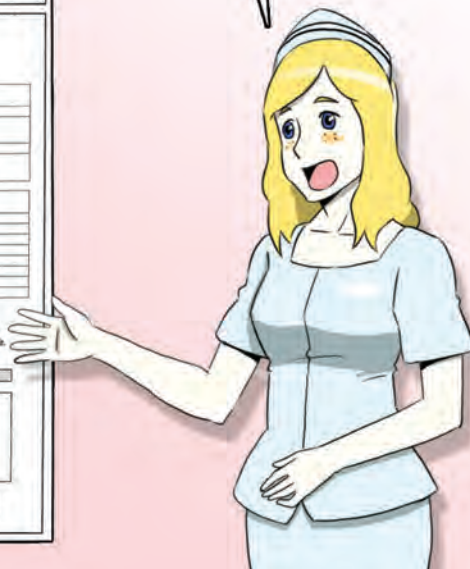

When using massage oil or water, or for a patient who is incontinent, we recommend using waterproof pads. These are available in 2 sizes : Professional(BMP) and mini.



Medical hyperthermia has effects on metabolism increase, good circulation, antibiotic treatment and especially detoxification.

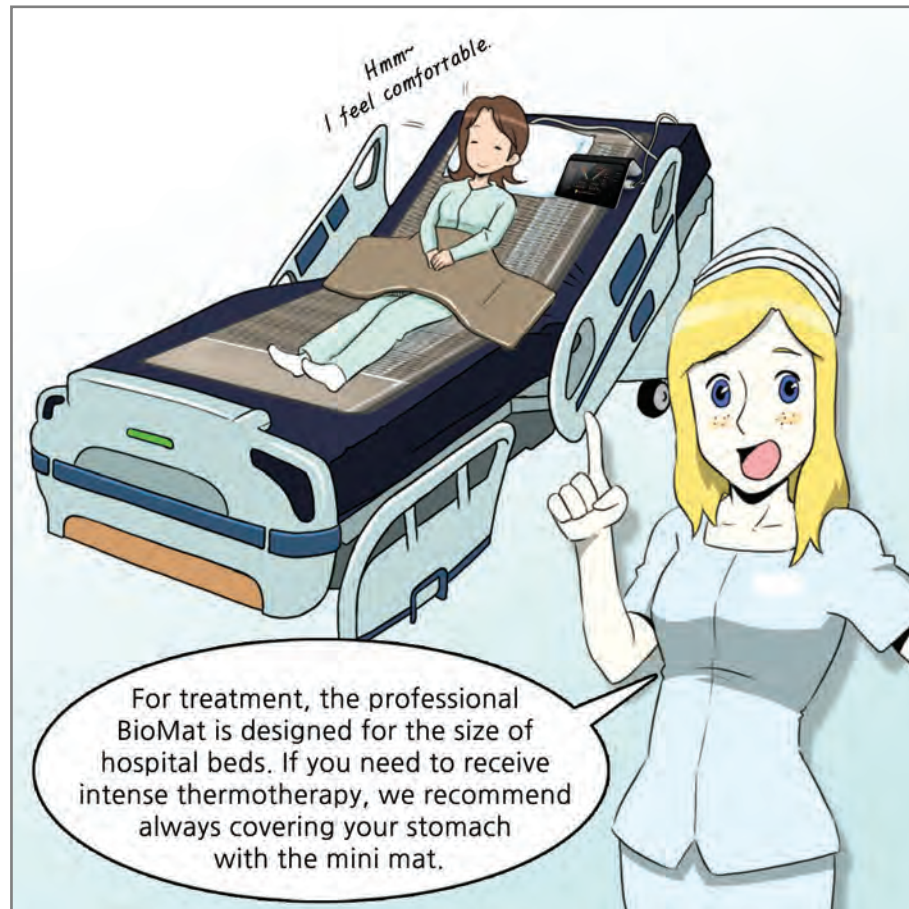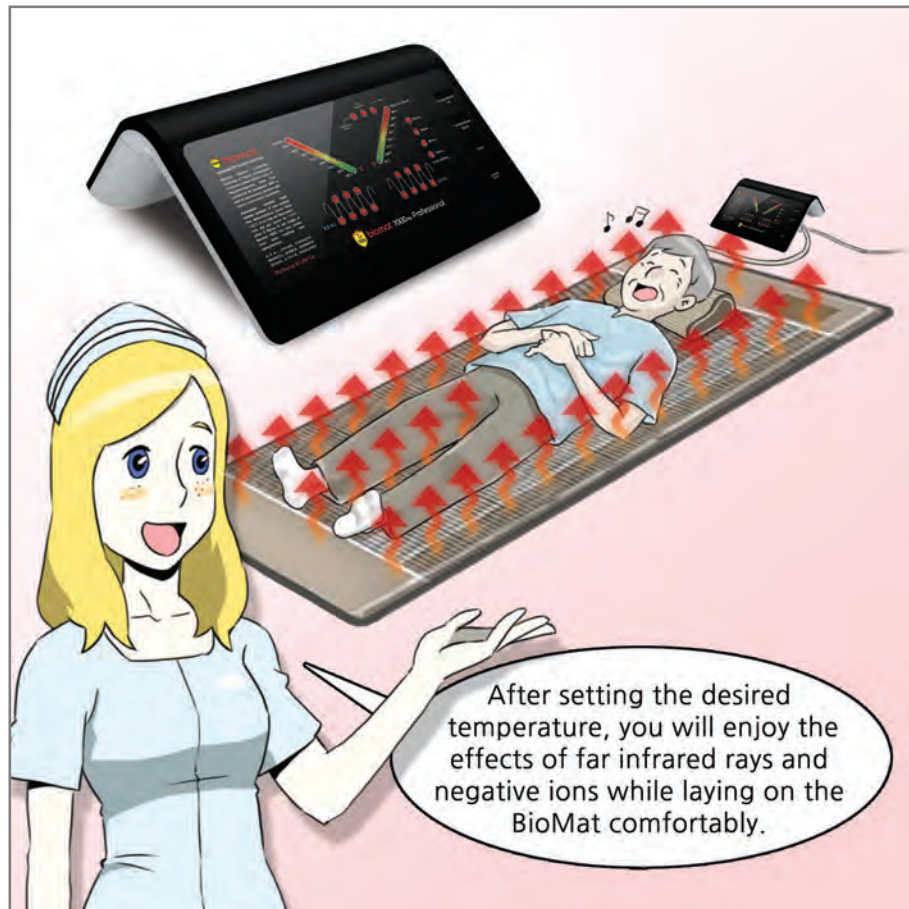

[About Chills No.1] Chills are the source of all kinds of disease.

How to use the mini-mat



The mini-mat is easy to carry.



[About Chills No.2] Chills on the body are worse in summer than winter.


You can use the mini-mat for animals.



Many people use the mini-mat for their pets. When using the mini-mat for a pet, we recommend using the waterproof pad.



Amethyst and tourmaline pillow emits natural far infrared rays.

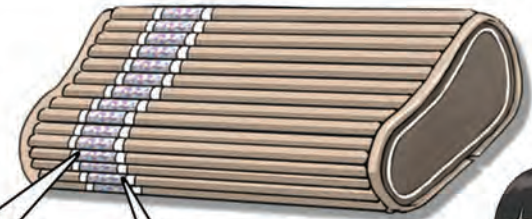

Using a Bio pillow filled with amethyst and tourmaline with the BioMat supplies natural far infrared rays and negative ions. Cranial benefits are many, including an enjoyable and comfortable deep sleep.

[About Chills No.3] Chilling the body makes it stiffer.

Prior to a high heat treatment,
drink a lot of water or Detoxi salt water.



Before your
thermotherapeutic treatment,
drink a lot of water or
Detoxi salt water to prevent
dehydration, as you may
sweat profusely.

It is important to consume Detoxi ionized salt
before performing thermotherapy.

IONIZED
SALT




Aids cells in controlling osmotic pressure
Strengthens kidney activity
Increases ionic activity inside the body

Detoxi 300 is ionized salt made by being heated for over 300 hours at 1,000°C. It is recommended that users take ionized salt with a high percentage of hydrochloric acid before using the BioMat because the salt can raise the therapeutic action of far-infrared radiation and help cells and kidneys to control osmotic pressure and to facilitate detoxification respectively.

[About Chills No.4] Chilling the body causes pain. The solution for the problem is only to warm the body.

How to use the BioMat for intensive treatment

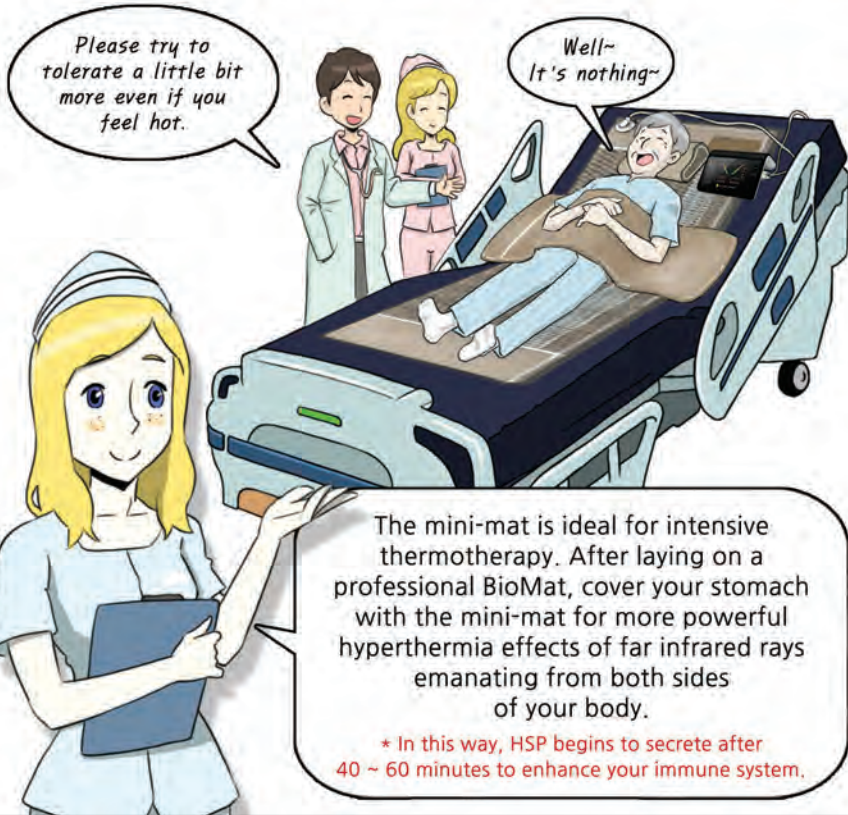


Please try to tolerate a little bit more even if you feel hot.

Well~ It's nothing~

The mini-mat is ideal for intensive thermotherapy. After laying on a professional BioMat, cover your stomach with the mini-mat for more powerful hyperthermia effects of far infrared rays emanating from both sides of your body.

* In this way, HSP begins to secrete after 40 ~ 60 minutes to enhance your immune system.

Generation of HSP

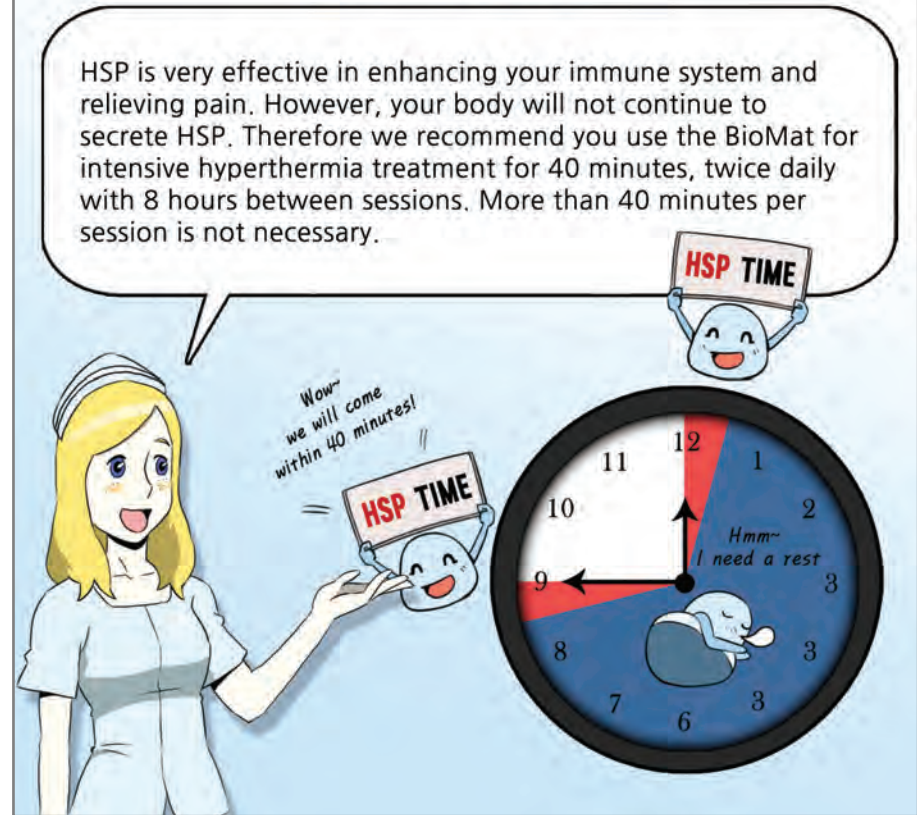

HSP is very effective in enhancing your immune system and relieving pain. However, your body will not continue to secrete HSP. Therefore we recommend you use the BioMat for intensive hyperthermia treatment for 40 minutes, twice daily with 8 hours between sessions. More than 40 minutes per session is not necessary.

Wow~ we will come within 40 minutes!

HSP TIME

HSP TIME


Hmm~ I need a rest

[About Chills No.5] Chilling the body causes obesity. The primary cause is that fat on the body becomes cold and stiffens.

When using the BioMat for intensive thermotherapeutic treatment vs. using it at low temperatures.


For an intensive hyperthermia session, lie on a professional mat and cover your stomach with a mini-mat. Set temperature at 70 °C, and stay in that position for 40 minutes.

In this way, HSP will be secreted. To sleep, set the temperature at 40 °C. If you use the BioMat while sleeping, your self-healing power will be enhanced more than 60 times.

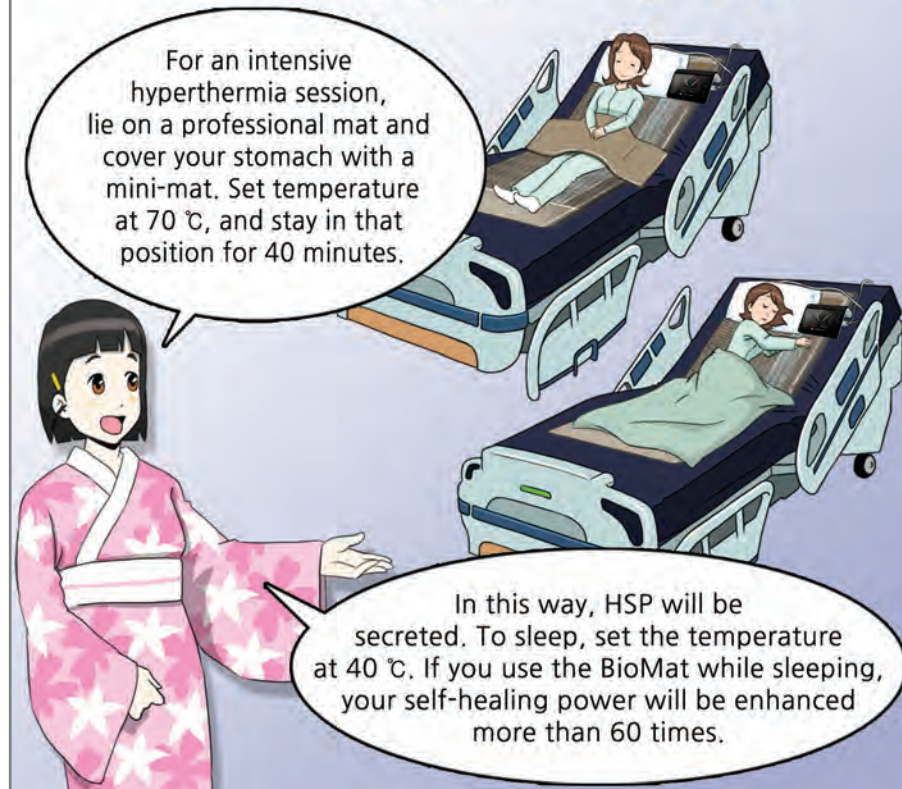

Precautions when using the BioMat.

For children and adults with diabetes or paralysis, pregnant women, breast-feeding mothers, people who have had bypass surgery or are using a pacemaker, always keep the BioMat setting at Low. If you have one of the conditions listed above, please consult your physician before using.

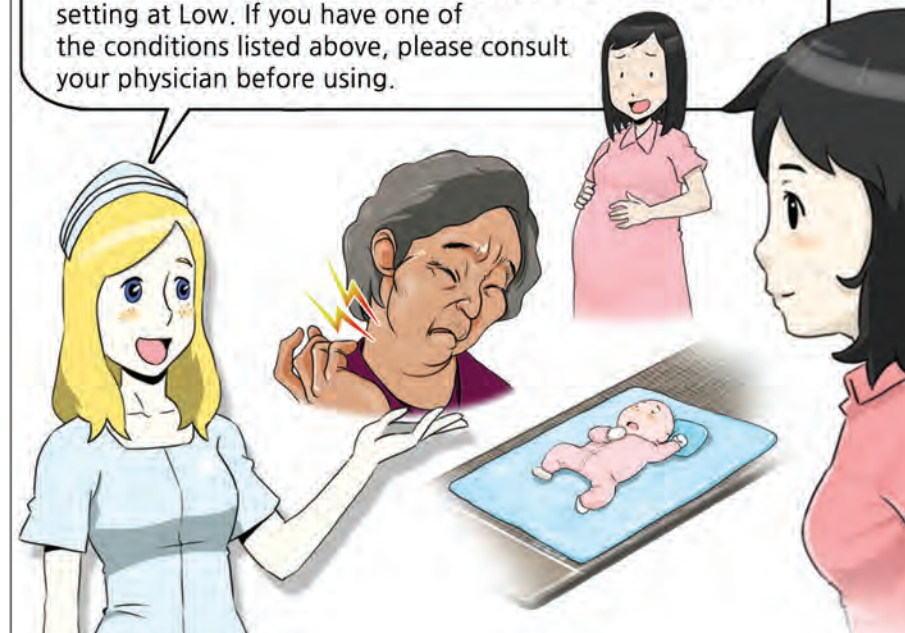

Children should be brought up warmly. (German Saying)

About the healing reaction



Toxic substances

I really hate far-infrared radiation ~

I want to live in a world without the Bio Mat ~

Me, too. Me, too!

Amethyst far-infrared radiation can have a radical effect on body temperature which can rise or fall by itself.

You may feel a lot of pain, or your physical condition may decline when toxins are excreted from your body, then your body will then begin to recover. That process is called the healing reaction. However, if you have a poor healing reaction, you can use the BioMat at a low temperature, or turn it off for a while and use it again later.

How to store the BioMat



Now ~


Let's see how to store the BioMat. Store your BioMat in a dry location to prevent the risk of electrocution due to damage from moisture.

The feeling of love begins with warmth. (Francois Sagan)

How to store the BioMat when not in use



Store the BioMat by folding it in the same way.

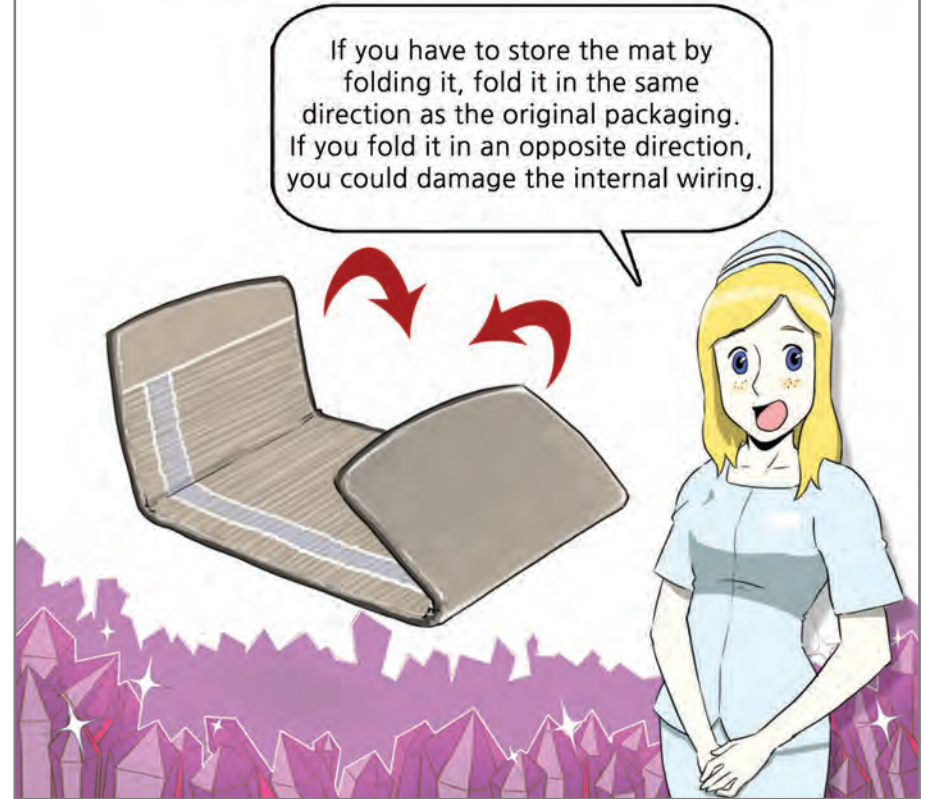

Low body temperature causes bad circulation due to the stiffening of red blood cells. This disturbs the process of excreting waste and toxic matters from the body and therefore causes various kinds of disease. (Naohisa Ishibashi, MD at the Kobe medical institute)

Do not store the BioMat by standing it upright.




You shall not wash or dry clean the BioMat.



Hyperthermia is effective if used for at least 30 minutes.

You can wash your quilted pad.

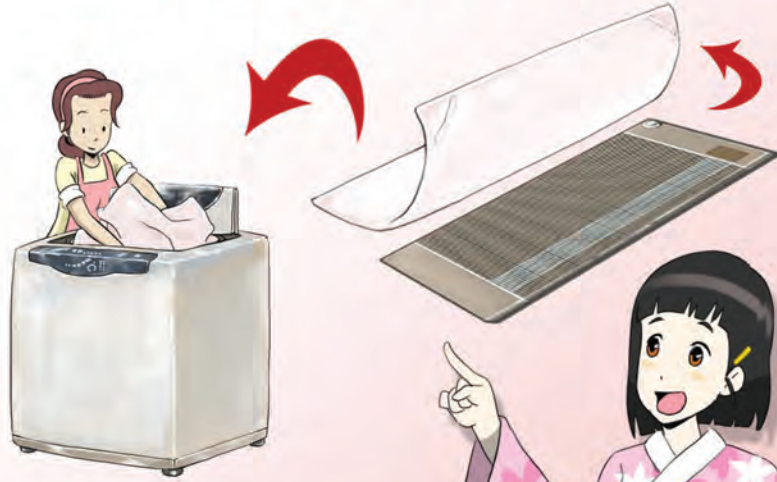

You can take the pad off from the BioMat and wash it separately, then hang the pad to dry.

Be careful when cleaning the controller.



When cleaning the controller, do not use wet towels. Water may cause a fire or damage the unit. Also do not use thinner, benzene, or wax. It may cause a change in shape or color.

Hyperthermia is a powerful weapon that the body uses against diseases.

You should not use any controller made by another company.

A controller from another company


To avoid risk of fire or other problems, use only the dedicated controller included in the package.

Do not place controller on top of BioMat or under covers.

We do not recommend placing the BioMat controller on the mat or in the bedding.

The metabolic rate increases by 100 percent whenever the body temperature rises by 1 degree. The increased metabolic rate elevates the production of interleukin, an anti-inflammatory substance. (Gaten, Doctor of physiology, the U.S.)

Do not jump on the mat.

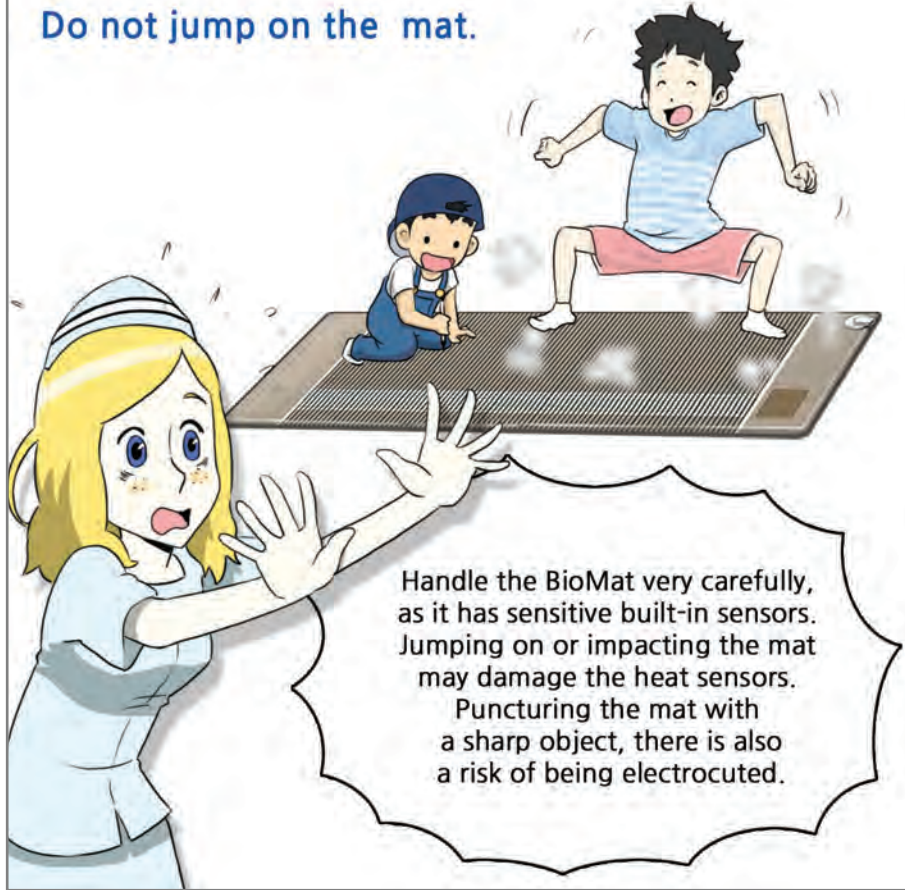

Handle the BioMat very carefully, as it has sensitive built-in sensors. Jumping on or impacting the mat may damage the heat sensors. Puncturing the mat with a sharp object, there is also a risk of being electrocuted.


Do not disassemble the unit.



If there is any problem with the product, do not try to disassemble or repair it. It may cause you to be electrocuted or other issues. If there is a problem, call Richway/The Life Corporation at (808)589-2800.

Hyperthermia, which cures diseases by raising body heat, allowed Dr. Cooley to win the Nobel Medical Prize in 1920.

Be careful not to entangle electric cord.

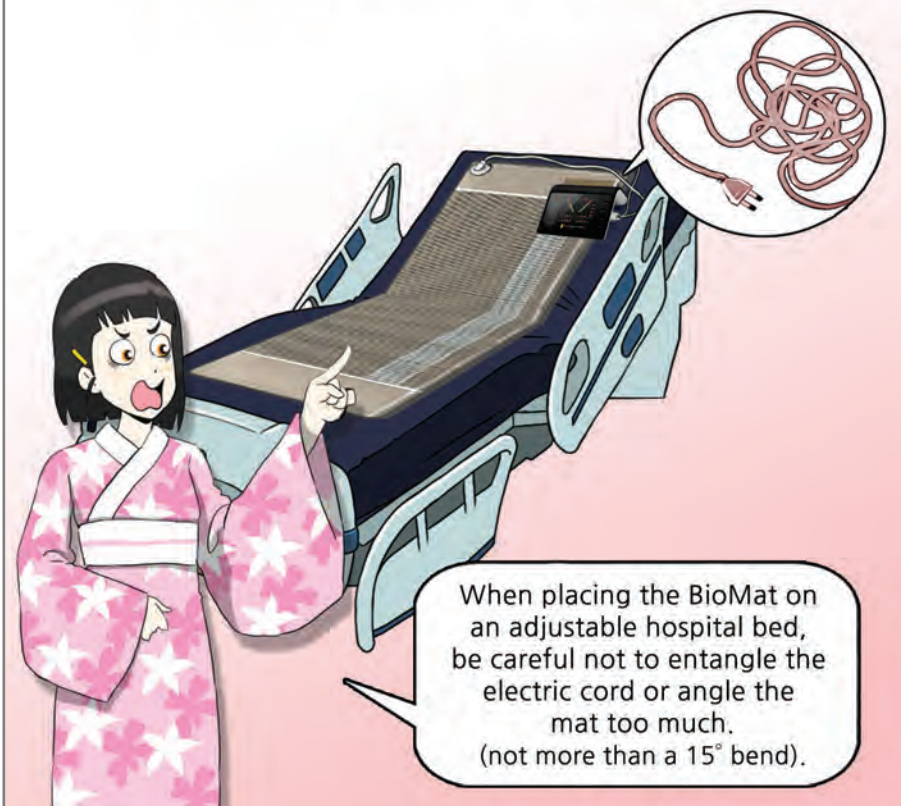


Do not put heavy objects on the BioMat.




Hyperthermia is the source from which a cure to all kinds of diseases can be found. (Parmenides, Doctor in Greece, 500BC)

Make certain your electrical cord is appropriate for your country's voltage

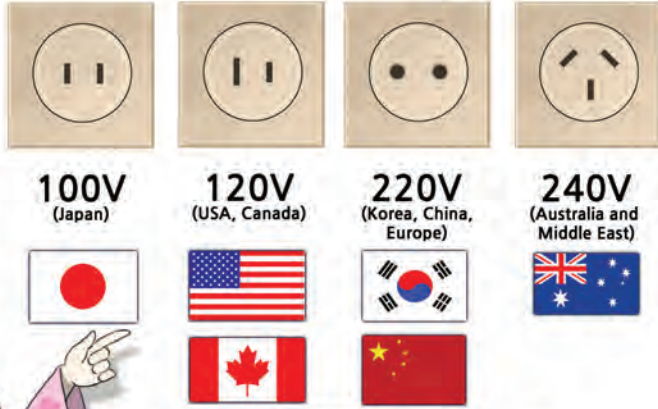

This product is manufactured at 240V, 220V, 120V, and 100V to meet the electrical standards of various countries. If you try to connect the product to a receptacle with different voltage, you may damage the product.

How to connect the power plug

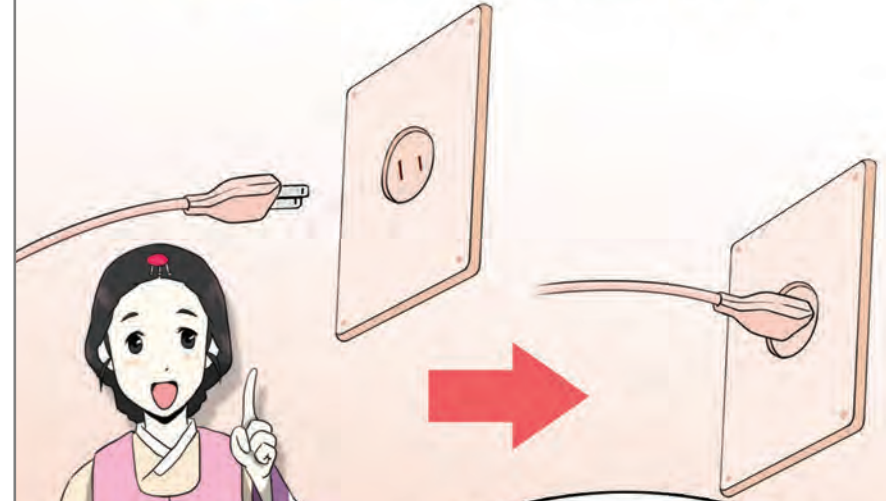


When connecting a power plug, you have to insert the plug tightly into the receptacle. If the plug is loose, the BioMat may not work properly due to lost contact.

Parmenides can be used safely to treat many kinds of diseases as a subsidiary therapy because it strengthens the body temperature and heartbeat. (Dr. Douglas Louis in the U.S., 1970)

Items that require attention when connecting a power plug



Do not insert the power plug with wet hands.

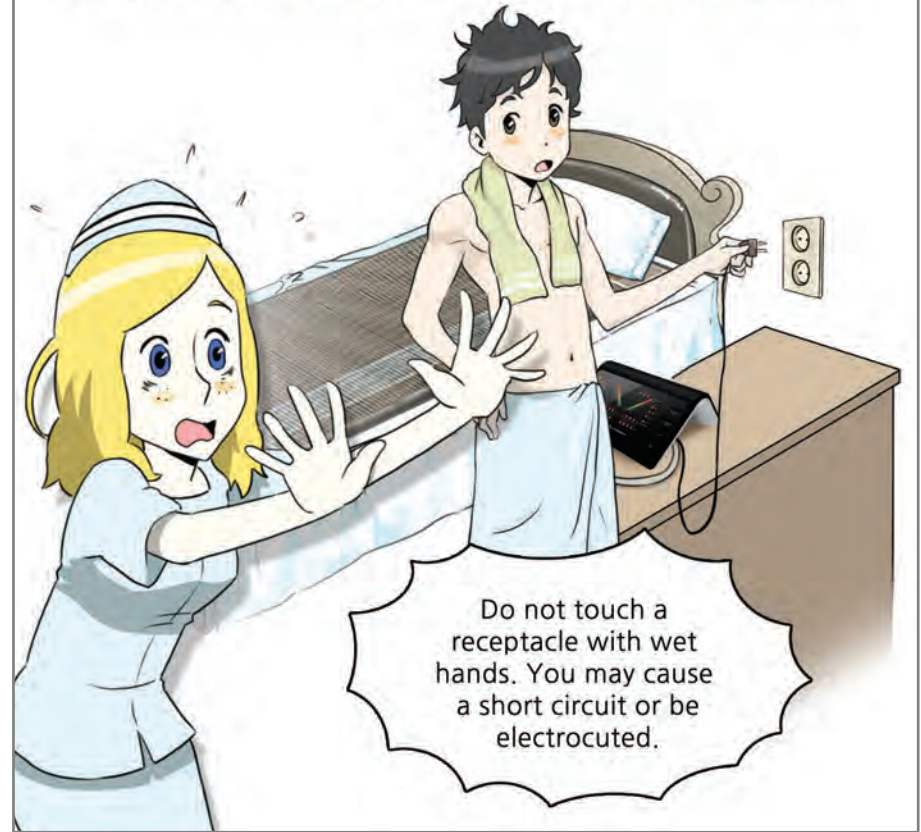

Hyperthermia is a kind of detoxification. It contains neutralizing pesticides, herbicides, and drugs, and it is an effective way to eliminate toxic matter from the system, initially through the skin, and secondarily through the intestines and the kidneys. (Dr. Gard, Cancer Support Center of the USA)



The BioMat awarded an acknowledgement from the United Nations.


The Richway BioMat received an award from the United Nations for its contribution to health care enhancement of the people around the world.



Richway received Pope Benedict XVI's Blessing of Richway for our contribution to healthcare products.

The Holy Father Benedict XVI cordially imparts the requested Apostolic Blessing to Dr. Jeanna Chun Kim, Richway International, Inc., for your receiving recognition by the United Nations for your leadership in promoting Health Care Enhancement and invokes an abundance of heavenly grace and the continued protection of the Blessed Virgin Mary.

On May 6, 2009. From Vatican

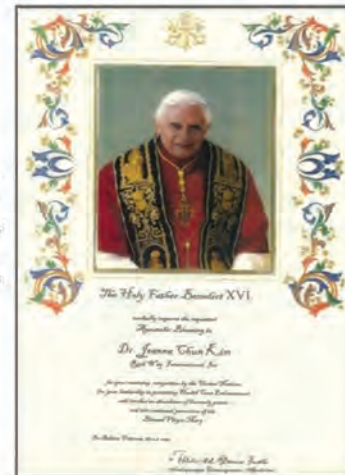

Richway received the blessing of Pope Benedict XVI for its contribution to save numerous HIV-infected children in Africa, and for treating many hard-to-cure diseases around the world.


We should call hyperthermia immunotherapy. (Dr. Denbayashi, GyohokuHospital)

Received recognition from Prince Husimi of the Japan Culture Agency.



Richway Japan received award from Prince Husimi of the Japan Culture Agency for its contribution to promote healthcare of the Japanese people.

Selected as a businessman of the year by the Wall Street Journal.

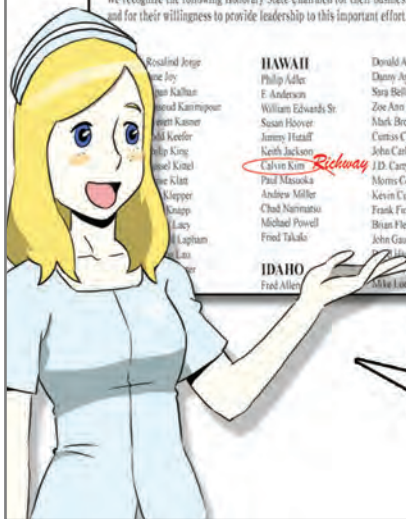

In 2006, Calvin Kim of Richway was selected as a businessman of the year by the Wall Street Journal, a prominent newspaper of the US.

It was confirmed that hyperthermia amplifies the effects of radiotherapy. (Ellen Jones, MD at Duke University Medical Center, The Wall Street Journal, 28 October 2007)

Richway was invited by the U.S. White House





To Chun and Calvin Kim
With Best Wishes,

Richway

In 2006, Richway was recognized as an excellent company selected by the President of the US, and invited to the White House for a dinner party hosted by President Bush.



The BioMat received the 510(k) approval for medical devices from the U.S. Food and Drug Administration

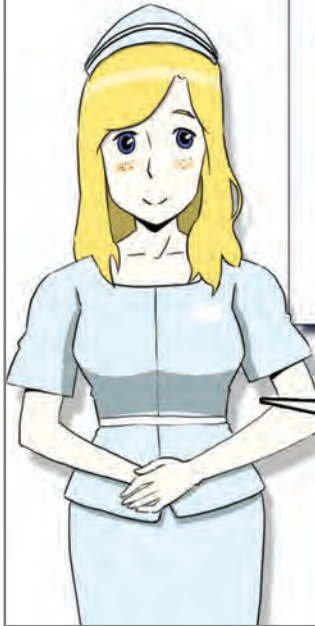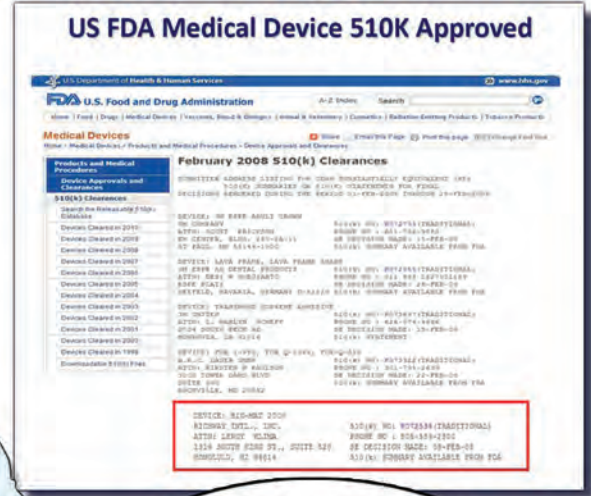

For the first time for a thermotherapeutic mat using far infrared rays, the BioMat received the 510(k) approval as a medical device from the US FDA. This signifies that the US FDA recognizes the safety and clinical effects of the BioMat.

In 1991, sixty percent of six hundred patients in Tokyo Women's Medical University got the positive effect of prolonging their life through 1000 instances of far-infrared radiation hyperthermia(treatment). (Yoo Yong-woon, the Head of Korea Far-infrared Radiation Medical Institute)

Attained certification from the Korea Food & Drug Administration, the Korea GMP, and the Korea Electrical Safety Corporation

The BioMat has been granted a medical device certification from the Korea Food & Drug Administration, and safety certification from the Korea Electrical Safety Corporation, and also received GMP certification, CE certification, and ISO 13485 medical device certification.



The BioMat : Approved by the Ministry of Welfare in Japan and the Japan Electrical Testing Laboratory




The BioMat has been approved by the Japan Ministry of Welfare, which applies the strictest criteria in the world as a medical device. It has also received a certification of safety from the Japan Electrical Testing Laboratory.




Hyperthermia can help to conquer any kind of disease. Diseases which cannot be treated with hyperthermia are irremediable. (Hippocrates, 450 BC)

The BioMat received a grand prize in the medical device field from the Royal Swedish Academy of Science which has appointed Nobel Prize nominees

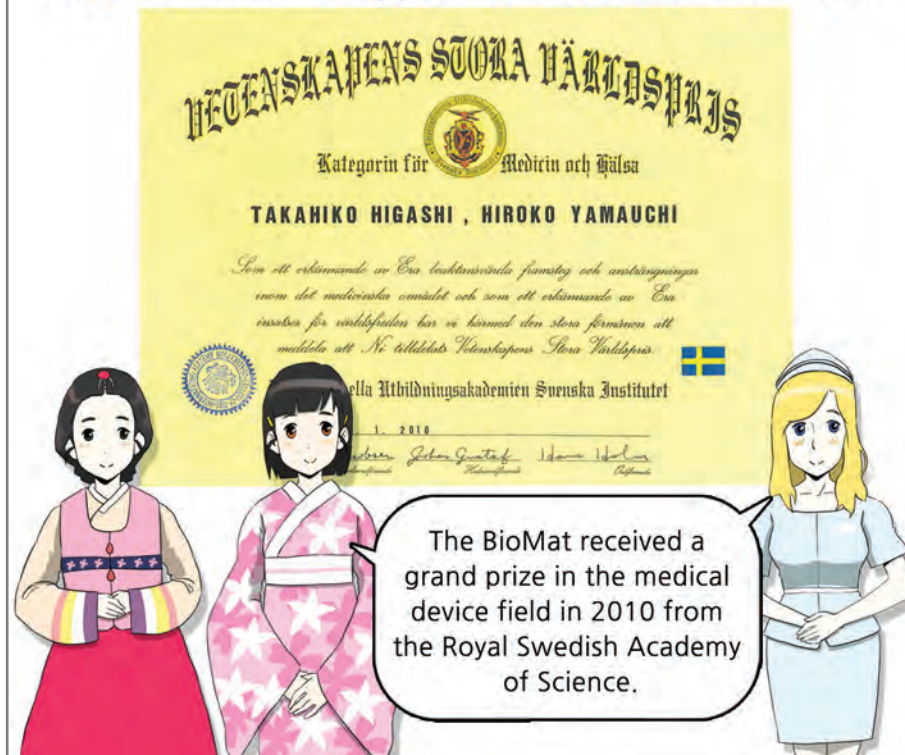

With the BioMat recognized by the world...

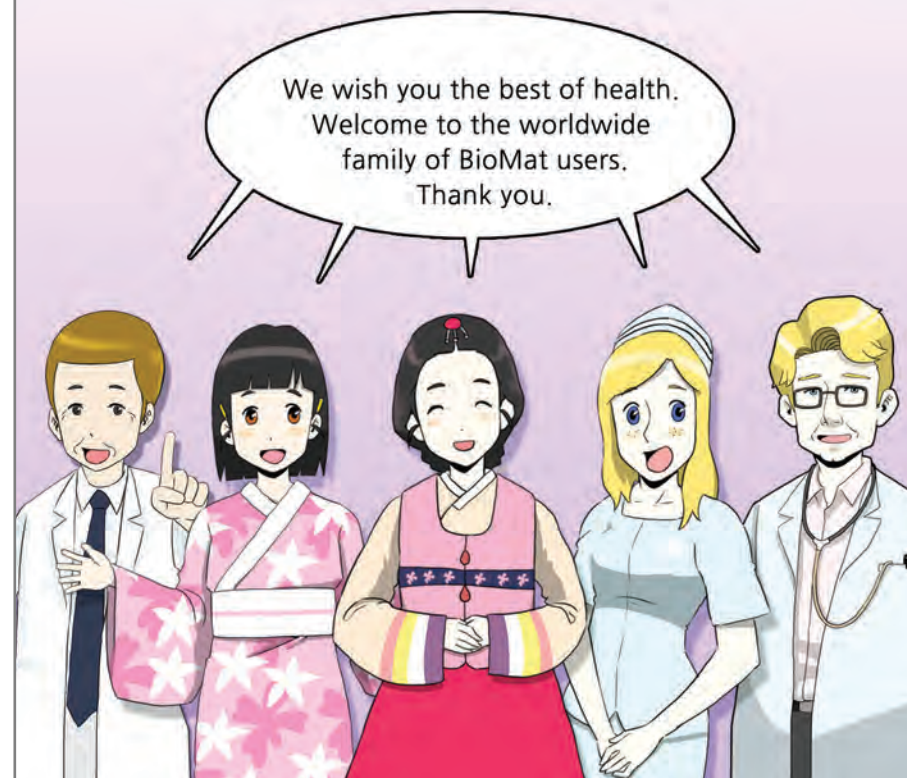

I will shine the disease-curing light upon you, and you will go outside and jump up and down like calves released from a stable. (Malachi 4:2)

You can use the BioMat when lying on your back or face down.

That looks comfy...

me too-me too-

I' m so jealous..

The BioMat is best used when you lie on your back~!

Don't you know lying face down on the BioMat is the best?

RICHWAY INTERNATIONAL / THE LIFE CORPORATION
LIFE WITH A BIOMAT

Happy 14th Birthday BioMat!

BIO-MAT

Happy Birthday CM

The power of recognition, activity, and proliferation, the three major parts of immunocyte, depend on maintaining a warm body temperature. Low body temperature causes the functions of immunocyte to cease and brings about a lethargic condition. (Ishihara, MD, 'The warm body saves your health')

R I C H W A Y

SPREADING THE MESSAGE AROUND THE WORLD



We believe that people deserve to be cured of their wounds,
renewed in spirit, spared from illness, enlightened from ignorance,
and saved from suffering. No one should be left out.
Our beauty does not lie in our appearance,
but in the passion of our showing and giving care.
Richway, the name you can trust....

